


HOW TO SECURE PERIPHERAL IV LINES



Step 1: Clean surface with disinfectant.
Allow to dry before using.

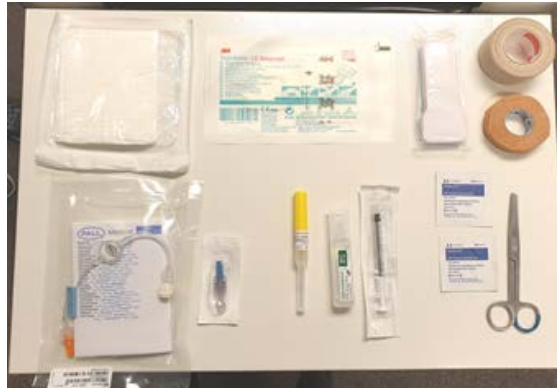


Step 2: Gather the equipment needed.



Step 3: Always use sterile drape underneath open equipment.



Step 4: Clean scissors with alcohol wipe.
Allow to dry before using.

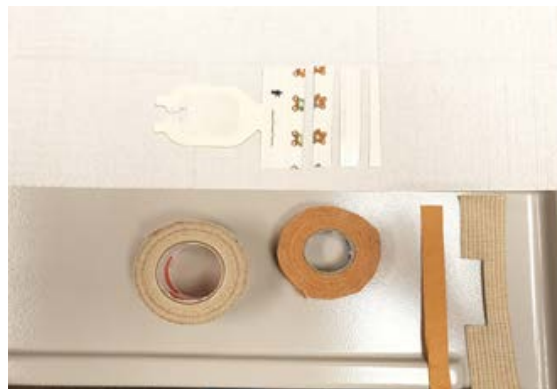

Step 5: Cut the steristrips then brown tapes



Step 6: Secure the IV leir with 2 x plain steristrips



Step 7: Cover with tegaderm. Ensure that the yellow hub is partially covered as this



Step 8: Attach 1 x teddy bear steri-strip under the hub.



Step 9: Attach another teddy bear steri-strip over the yellow hub.



Step 10: Attach steri-strip with written date to the bottom of the tegaderm. Connect the extension. Attach the Elastoplast at the bottom end of the board.



Step 12: Attach a long strip of brown tape under the hub and around the fingers



Step 13: Reinforce with another long strip of brown tape above the hub.